


Rocket® FBS Sampling Wand



The **Rocket® FBS Sampling Wand** combines both a fetal scalp blade and a capillary sample tube in one simple instrument.

It's as easy to use as a ball point pen!

- **IMPROVES EASE OF SAMPLING** by avoiding the need to change between separate FBS blade and capillary tube holders whilst visualising the fetal scalp.
- **SIMPLE, SAFE BLADE LOADING AND RETRACTION:** The pre-loaded FBS Blade is advanced for use by depressing the rear moulding after the scalp incision. It is automatically retracted by depressing the button on the wand shaft, once retracted the capillary tube is ready for sampling.
- **RELIABLE SAMPLING:** In conjunction with the world's leading manufacturers of blood gas analysers, Rocket Medical have specifically developed Sampling Wands and complete FBS kits compatible with the following systems:
 - ✓ Radiometer® ABL™ 70, 90, 800, 700, 77, and 5/500 Series
 - ✓ Siemens® (Bayer) RapidLab™ Series
 - ✓ Instrumentation Laboratory® GEM™
 - ✓ ROCHE micro volume systems.
- **COMPLETE WITH SAMPLING TUBES:** Each Sampling Wand comes complete with an FBS Blade and five dedicated capillary tubes in one of four volumes.
- **STERILE, FOR SINGLE PATIENT USE,** supplied in cartons of 10 wands.

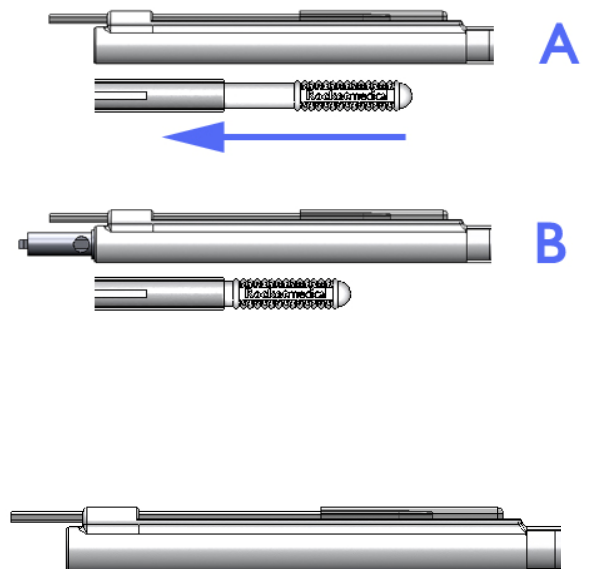

Description	FBS Kit equivalent	Volume	Code
Rocket® FBS Sampling Wand complete with 5 heparinised capillary sampling tubes, sterile for single patient use, supplied in cartons 10 wands.	R57018	55µl	R57021-55-SW
	R57012 R57013	75µl	R57021-75-SW
	R57024	85µl	R57021-85-SW
	R57026	100µl	R57021-10-SW

