



Rocket® BulbTip™ & BulbTip Echo™ Embryo Transfer Catheters

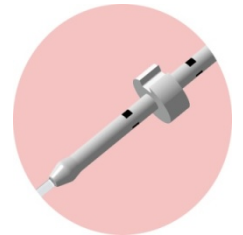

Rocket® BulbTip™ & BulbTip Echo™ catheters are the latest generation of embryo transfer products from Rocket Medical. Designed with new features to further improve ease of insertion and minimize uterine trauma, they continue our long established reputation for continuous product development.

- **NEW ECHO BANDED VERSION** to improve visualisation under ultrasound guidance.
- **ULTRAFINE INNER TRANSFER CATHETER** to minimise media volume on transfer. Externally supported to ease loading and to aid insertion into the sheath during transfer. Distal 60mm unsupported.
- **BULB TIP OUTER SHEATH** with smooth bulb formed tip to improve location during insertion, reduce cervical trauma and to guide the ultra fine inner catheter correctly into the uterine cavity. Particularly suitable for patients where cervical access is difficult.
- **SLIDING SILICONE MARKER** to allow insertion depth to be pre-set prior to transfer.
- **PATENTED FINGER GRIP** gives unequalled grip and control while ensuring accurate placing of the catheter tip into the uterine cavity.

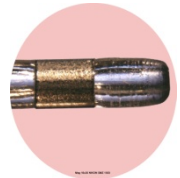

Embryon® BulbTip™ Range	Length	Code
Rocket® BulbTip™ Embryo Transfer Catheter Sterile, for single patient use, supplied in cartons of 10 units	18cm	R57635-00-18
	23cm	R57635-00-23
Rocket® BulbTip Echo™ Embryo Transfer Catheter Sterile, for single patient use, supplied in cartons of 10 units	18cm	R57635-EC-18
Rocket® BulbTip™ Trial Transfer Catheter Sterile, for single patient use, supplied in cartons of 10 units	18cm	R57636-00-18
	23cm	R57636-00-23
Rocket® BulbTip™ Stylet Sterile, for single patient use, supplied in cartons of 10 units	18cm	R57591-BT-18
	23cm	R57591-BT-23

